


INTENDED USE: For the collection and disposal of aspirated fluids.

DIRECTIONS FOR USE

1. Place lid on canister and press firmly around entire perimeter.
2. Apply pour spout cap firmly over pour spout.
3. Firmly attach vacuum tube to vacuum port and attach patient tube to patient port.
4. Check all connections for proper seal. Test the assembly for vacuum leaks by turning on vacuum source and occluding the patient tubing with finger or thumb.
5. To prevent fluid contamination to vacuum system, make sure vacuum line is attached to vacuum port.

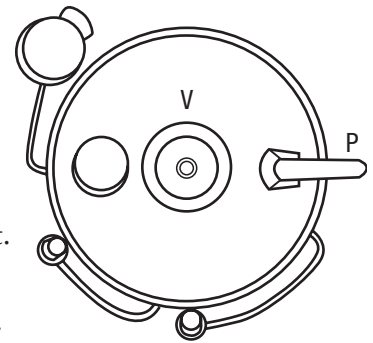


P = PATIENT PORT
V = VACUUM PORT

DISPOSAL

1. Turn off vacuum source and disconnect all tubing.
2. Seal vacuum and patient ports with attached port caps.
3. Remove canister from bracket and transport to disposal area. Do not lift canister by lid. The weight of the contents may cause the lid to separate from the canister.
4. Dispose of according to hospital policy.

⚠ CAUTION

- **Single Use Only. Do not attempt to clean, sterilize or reuse canister.**
- Do not exceed vacuum level of 25" (640mm) Hg.
- Not recommended for use in liposuction procedures.
- Do not apply continuous vacuum longer than 24 hours.
- Not intended as a measuring device, only for general reference, not specific measurement.
- Canister contents are considered potentially hazardous, use appropriate PPE and handle accordingly.
-  Check expiration date on canister for use-by date.
- Store in a dark place. Long term exposure to light may compromise product performance and result in breakage during use.

Loss of Vacuum: Check that vacuum is on, canister is properly sealed and that all connections are tight and tubing is not kinked. If loss continues, replace canister.

Report any product malfunctions or complaints to Bemis Manufacturing Co. using the contact information below.

See brochure for complete line of canisters, hardware and accessories.



Bemis Manufacturing Co., 300 Mill Street, Sheboygan Falls, WI 53085 • 920-467-4621 • 800-558-7651 • FAX 920-467-8573
Bemis Limited, Farrington Road, Rossendale Road Industrial Estate, Burnley Lancashire BB1 5SW, England • 44-1282-438717 • FAX 44-1282-412717
e-mail HCG@BemisMfg.com • www.bemishealthcare.com

