


QUICK-CONNECT II REGULATOR BRACKET

DESIGNED FOR FAST, SAFE AND SECURE SUCTION CANISTER CONNECTIONS



EFFICIENCY

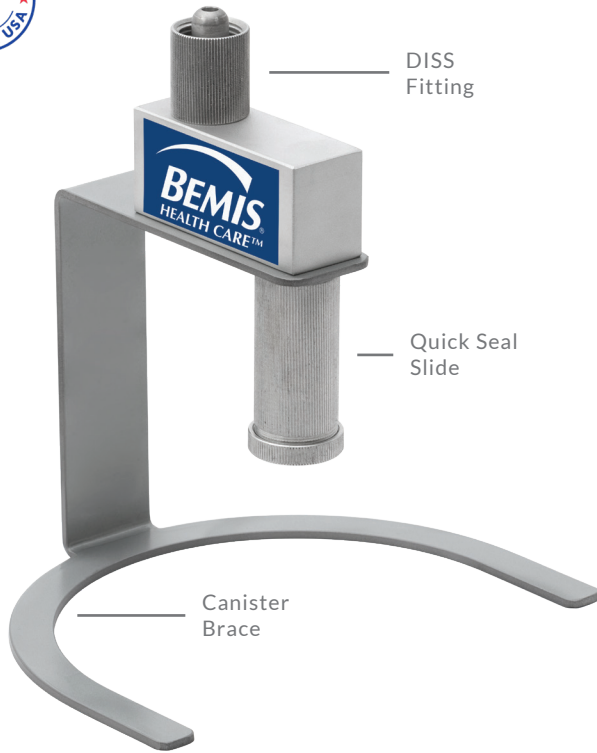

SAFETY



DURABILITY

FEATURES + BENEFITS

The Quick-Connect II is a direct-to-vacuum regulator bracket designed for fast, safe and secure suction canister connections — every time. For use in any patient care area with a central vacuum wall regulator.

EFFICIENCY

- **Space-saving design** — Clinician-inspired solution hangs in-line from regulator to reduce wall clutter
- **Simple installation** — DISS connection means no tools or drilling needed; eliminates need for wall plates, regulator connection tubing and wall mount accessories
- **Intuitive Quick Seal Slide** — Provides streamlined ergonomics for fast canister setup and use

SAFETY

- **Fail-safe design** — Ensures proper vacuum hook-up for reliable suction; requires minimal training for safe use
- **Secure canister** — Unique Quick Seal Slide locks canister onto brace and assures airtight seal
- **Safe port capping** — Canister Brace allows stable port capping for safe disposal

DURABILITY

- **Reusable bracket** — Made with rugged, quality materials for repeated use
- **Minimal maintenance** — Easy, in place cleanability and disinfecting

QUICK-CONNECT II REGULATOR - BRACKET CONFIGURATIONS

800CC CANISTER CONFIGURATIONS



800cc Hi-Flow SOLO
#534110 / 8002 055

1200CC CANISTER CONFIGURATIONS



1200cc Hi-Flow
#534310 / 484410

QUICK-CONNECT II CONFIGURATIONS + COMPATIBILITY

PRODUCT/DESCRIPTION	CATALOG#	PACKAGING
---------------------	----------	-----------

800CC CONFIGURATIONS



QUICK-CONNECT II 800CC SOLO REGULATOR BRACKET

534110

1/CASE

- For 800cc Hi-Flow SOLO canisters (#8002 055)

HxD without canister: 7.188in x 4.625in

HxD with canister: 11.875in x 5.75in



800CC HI-FLOW SOLO CANISTER

8002 055

70/CASE

- Float valve shutoff with Aerostat filter

Bottom Diameter: 3.50in Top Diameter: 5.375in Height: 7.25in

1200CC CONFIGURATIONS



QUICK-CONNECT II 1200CC REGULATOR BRACKET

534310

1/CASE

- For 1200cc Hi-Flow (#484410) or Sure-Flow canisters (#634410)

HxD without canister: 7.00in x 5.00in

HxD with canister: 12.75in x 6.50in



1200CC HI-FLOW CANISTER

484410

48/CASE

- Float valve shutoff with Aerostat filter, two elbows

Bottom Diameter: 4.125in Top Diameter: 6.25in Height: 7.00in

Note: Vacuum regulator requires the use of a male N.P.T. X DISS fitting to use Quick-Connect II bracket. If your facility does not currently use DISS fittings on regulators, refer to Catalog# 531910. Bemis Health Care does not manufacture or distribute vacuum regulators. For installation, use, cleaning and maintenance instructions, refer to Quick-Connect II IFU (5990001).

Bemis Manufacturing Company

Health Care Division | 300 Mill Street | PO Box 901 | Sheboygan Falls, WI 53085 USA

Email: Customer.Service@BemisHealthCare.com | Phone: 920.467.4621 | Fax: 800.499.8160

4/2025 Rev 02

©2025 0B77522116

