

EXPIRATION INFORMATION: CANISTERS & LINERS

EXPIRATION DATE INFORMATION

- All Bemis Canisters and Liners include an expiration date on the carton label and individual product.
- Expiration is three (3) years from the date of manufacture.
- Canisters and Liners should not be used past the labeled expiration date.
 - Structural integrity can be compromised by factors such as, heat, time, and light exposure.
 - Placing a negative vacuum on a compromised canister/liner can result in product failure.

BEMIS RIGID CANISTERS

Expiration date can be found on:

- The outer product carton label
- Individual canisters: notated on front of canister by: ISO symbol | year/month/day
- Please note there is a manufacturing date printed on the bottom side of all Contour and Hi-Flow lids, not to be confused with the expiration date

RIGID CANISTER PRODUCTS

- **Hi-Flow™** (800cc, 1200cc, 2000cc, 3000cc)
- **Hydrophobic** (800cc, 1200cc, 2000cc)
- **Contour™** (2000cc, 3000cc)

ON-PRODUCT EXAMPLE



BEMIS LINERS


Expiration date can be found on:

- The outer product carton label
- Individual liners: stamped on lid top notated by: Catalog #, manufacturing date, expiration date

LINER PRODUCTS

- **Quick-Fit™** (1000cc, 1500cc, 3000cc)

ON-PRODUCT EXAMPLE

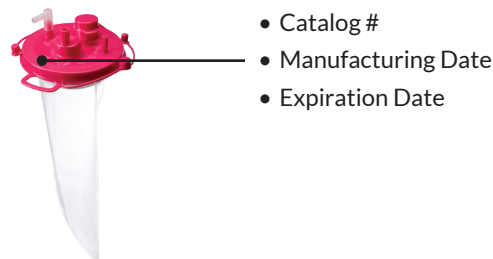

Bemis Manufacturing Company

Health Care Division | 300 Mill Street | P.O. Box 901 | Sheboygan Falls, Wisconsin 53085 USA

Email: Customer.Service@BemisHealthCare.com | Phone: 920.467.4621 | Toll-Free Phone: 800.558.7651 | Fax: 800.499.8160

6/2024 Rev 02

©2024 0B7752115

